FETAL ILCP BIOMARK

General informations

Samples:

2 x 48 well plates

24 single cells LTiP

24 single cells ILCP

22 single cells $\alpha_4\beta_7$

6 well of 10 cells of each

Genes:

66 Solaris Assay

27 Taqman Assay

93 Assay in total

Biomark runs:


1 96x96 chip (Solaris Assays)

2 48x48 chips (Tagman Assays)

Analysis:

Single cells considered positive if at least 3 housekeeping genes are expressed Each gene is analyzed by Δ Ct method (with average of 3 housekeeping genes Actb, Hprt and Gapdh)

Overall results

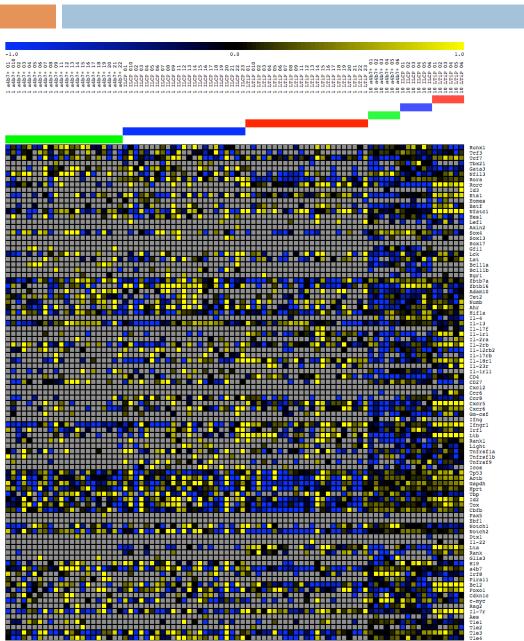

Analyzed single cells: 23/24 single cells LTiP (Red) 23/24 single cells ILCP (Blue) 22/22 single cells $\alpha_4\beta_7$ (Green) +6 wells of 10 cells of each

Legend:

Columns: Sample

Row: Assay

Grey: Not Detected

Blue: Low expression

Yellow: High Expression

Results are presented as $-\Delta Ct$, each Gene expression is centered around median and at similar scales.

